UTS CAD

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Seite 3

Operating Instructions

Page 15

Mode d'emploi

Page 27

Istruzioni d'uso

Pagina 39

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Página 51

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Página 63

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

Table of Contents

1 Introduction / Signs and Symbols	17
1.1 Preface	17
1.2 Introduction	17
1.3 Notes regarding the Operating Instructions	17
2 Safety First	17
2.1 Indications	17
2.2 Health and safety instructions	17
3 Product Description	18
3.1 Construction and functional description	18
3.2 Indications, Contraindications	18
4 Installation and Initial Start-Up	19
4.1 Unpacking and checking the contents	19
4.2 Installation and initial start-up	19
5 Handling, Operation	19
5.1 Fixture to determine the CE angle value	19
5.2 Fixture to determine the BP angle value	19
5.3 CE scale	20
5.4 BP scale	20
6 Practical Use on the Patient	20
6.1 Preparing the registration	20
6.2 Performing the registration	21
7 Maintenance, Cleaning, Diagnosis	23
7.1 Monitoring and maintenance	23
7.2 Cleaning	23
7.3 Maintenance	23
8 What If ...	24
8.1 Technical malfunctions	24
8.2 Repair	24
9 Product Specifications	24
9.1 Delivery form	24
9.2 Technical data	24
10 Miscellaneous	25

Views of UTS CAD / List of Parts

List of Parts

- 1 Basic bow
- 2 Bite fork
- 3 Bite fork & centric tray adapter
- 4 BP screw
- 5 BP indicator
- 6 BP scale
- 7 CE screw
- 8 CE indicator
- 9 CE scale
- 10 Screw bite fork fixation
- 11 Carrier
- 12 Rods
- 13 Clasp
- 14 Clasp washer
- 15 Clasp washer

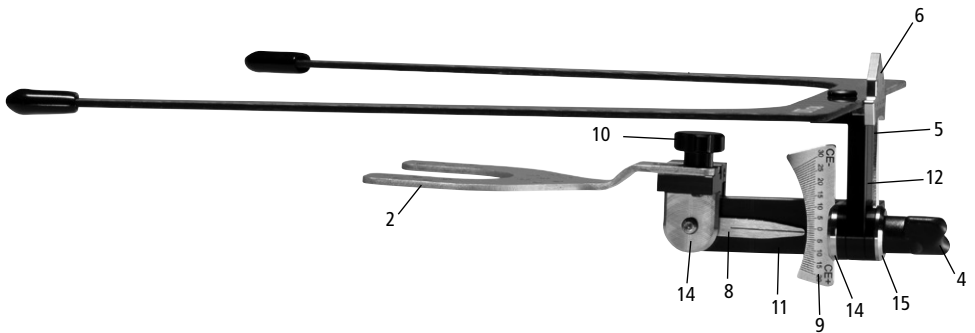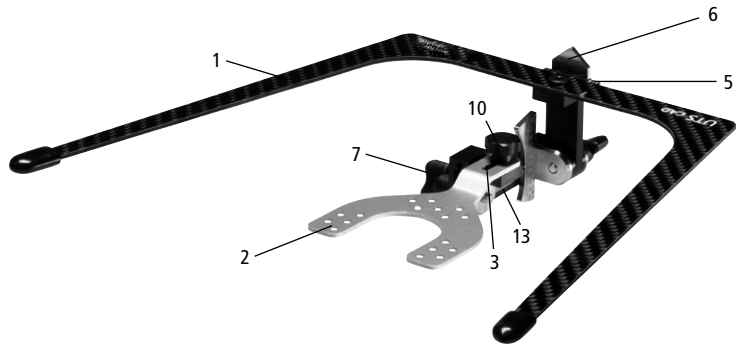1 Introduction/Signs and Symbols

2 Safety First

1.1 Preface

Dear Customer,

Thank you for having purchased the UTS CAD. It is a highly-technical quality instrument. This product is distinguished for its good quality and high precision.

The UTS CAD has been designed according to the latest industry standards. Inappropriate use may damage the equipment and be harmful to personnel. Please observe the relevant notes and read the Operating Instructions carefully.

Enjoy working with the UTS CAD.

1.2 Introduction

The signs and symbols in these Operating Instructions facilitate the finding of important points and have the following meanings:

Important information

Risks and dangers

Contraindication

1.3 Notes regarding the Operating Instructions

Furnace concerned: UTS CAD
Target group: Dentists, orthodontists, dental technologists

These Operating Instructions facilitate the correct, safe and economical use of the UTS CAD.

Should you lose the Operating Instructions, extra copies can be ordered from the respective Service Centre or downloaded as PDF from the internet.

We recommend this chapter to all individuals working with the device or perform cleaning and maintenance work.

2.1 Indications

When manipulating UTS CAD make sure not to injure the eyes of the patient with the components, particularly with the end of the bow.

Risks and dangers:

When manipulating UTS CAD make sure not to injure the eyes of the patient with the basic bow.

2.2 Health and safety instructions

UTS CAD must only be used for the purpose described in Chapter 3. Further instructions to assure proper use of UTS CAD:

- The instructions, regulations, and notes in these Operating Instructions must be observed.
- UTS CAD must be properly maintained (see Chapter 7).

3 Product Description

3.1 Construction and functional description

The UTS CAD is a registration device for the determination of the angles from the occlusal plane to Camper's plane and from the occlusal plane to the bipupillary line.

The measured angles and/or deviation of the occlusal plane from CE/BP can be transferred to the denture design software, which then reproduces the correct virtual position of the occlusal plane for the design of the denture.

CE = Camper's plane

Camper's plane runs from the nasal spine to the lower edge of the external auditory meatus.

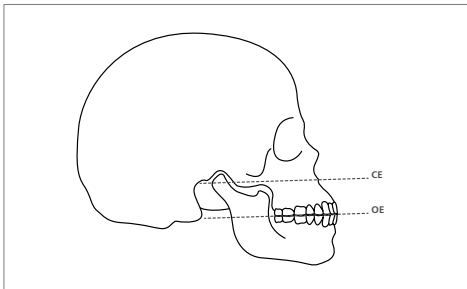
Image: OP to CE

BP = Bipupillary line

The bipupillary line runs through the centre of the pupils of the left and right eye of the patient.

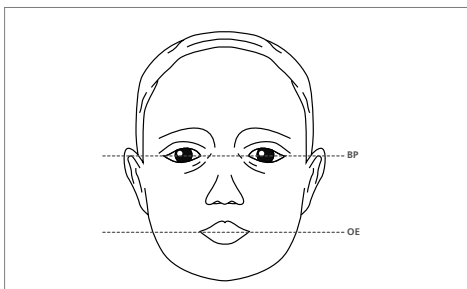
Image: OP to the bipupillary line

3.2 Indications, Contraindications

Indications:

Determining the angle of the occlusal plane to Camper's plane

Determining the angle of the occlusal plane to the bipupillary line

Contraindication:

The device cannot replace a transferbow for skull/TMJ-related model transfer to the articulator.

4 Installation and Initial Start-Up

4.1 Unpacking and checking the contents


Please make sure that the contents are complete.

- 1 UTS CAD basic bow, complete
- 1 3D bite fork
- 1 Operating Instructions


4.2 Assembly and initial start-up

The UTS CAD is supplied preassembled.

You only have to mount the bite fork (2) at the bite fork & centric tray adapter (3) and secure the bite fork fixation with the screw (10).

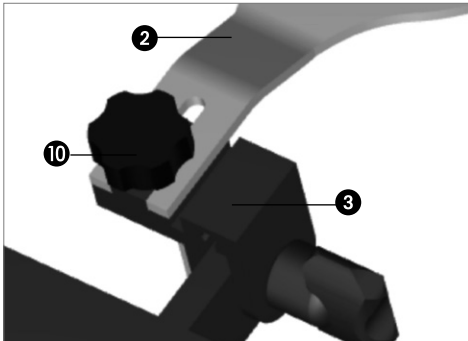
Please make sure that the bite fork (2) is correctly slid into the bite fork adapter (3) until it stops and that it is flush. Only then is it ensured that the bite fork & centric tray (2) can be safely and correctly secured at the bite fork & centric tray adapter (3).

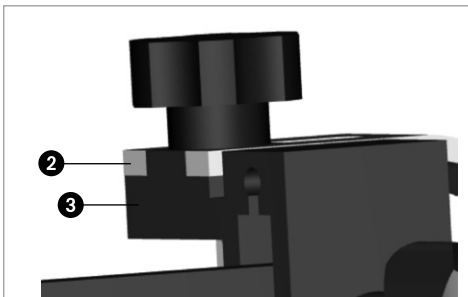5 Handling, Operation


5.1 Fixture to determine the CE angle value

The joint can be unlocked by loosening the CE screw (7) at the clasp (14). Then the joint can be freely swivelled. After adjustment, the joint can again be locked with the CE screw (7).

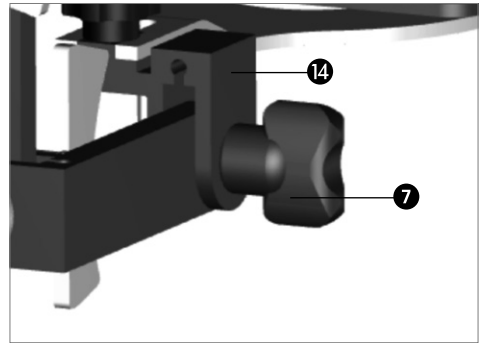5.2 Fixture to determine the BP angle value


The joint can be unlocked by loosening the BP screw (4) at the clasp (14). Then the joint can be freely swivelled. After adjustment, the joint can again be locked with the BP screw (4).

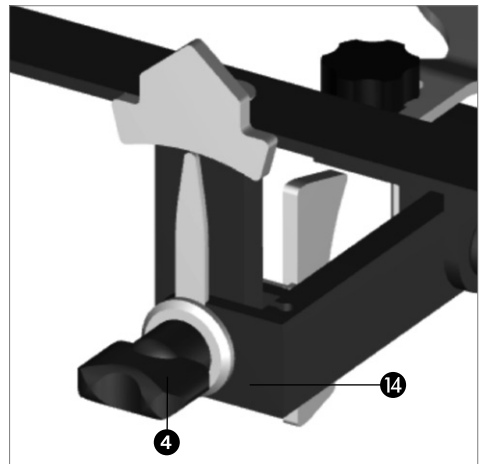6 Practical Use on the Patient


5.3 CE scale

The CE scale (9) permits indications from +20 to -30°. The CE scale (9) is divided into a plus and a minus side. The scale is readable with an accuracy of 1°. In this example the CE indicator (8) shows a CE angle of +2°.

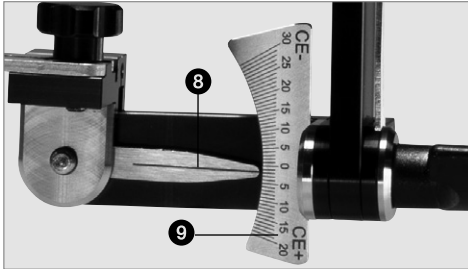
6.1 Preparing the registration

Place the bite fork (2) between the two Gnathometer CAD bite rim supports and instruct the patient to keep the mouth closed so that the bite fork is secured in a stable manner.

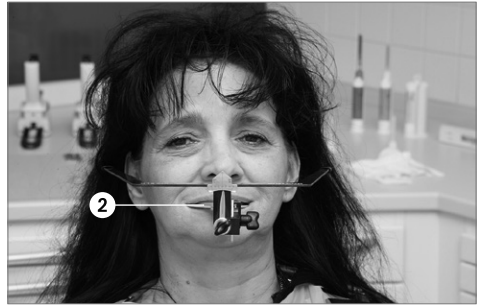5.4 BP scale


The BP scale (6) permits indications of +/- 15°. The BP scale (6) is divided into a plus and a minus side of 15°. The scale is readable with an accuracy of 1°. In this example the BP indicator (5) shows BP angle of 0°.

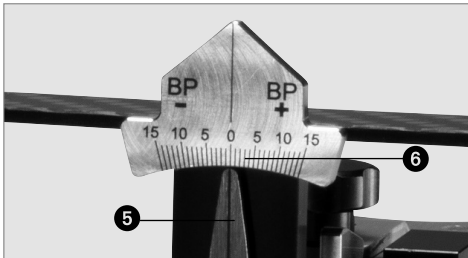6.2 Performing the registration


6.2.1 Mounting the transferbow

The patient should be in an upright position. Now, connect the UTS CAD with the bite fork (2) or the centric tray. For that purpose, slide the bite fork (2) or the centric tray onto the bite fork & centric tray adapter (3) and secure the components with the bite fork fixation screw (10).

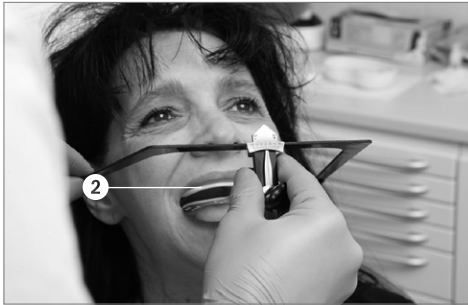
6.2.3 Aligning the transferbow with the BP line

The joint can be unlocked by loosening the BP screw (4) at the clasp (14). Now, align the front side of the basic bow (1) parallel to the BP plane. After alignment, the joint can again be locked with the BP screw (4).

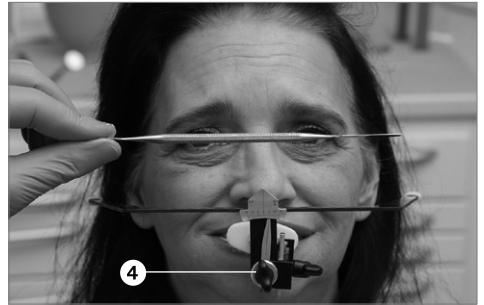6.2.2 Aligning the transferbow with CE

The joint can be unlocked by loosening the CE screw (7) at the clasp (14). Align the lateral arms of the basic bow (1) parallel to CP. If the line alignment is correct, tighten the CE screw (7). This position is now secured.

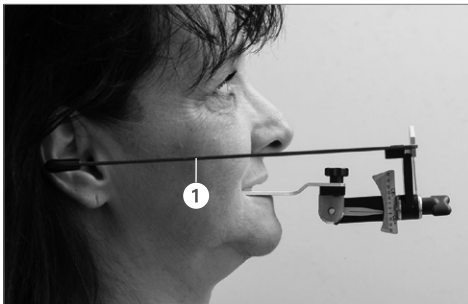
6.2.4 Removing UTS CAD from the patient

Instruct the patient to open the mouth. Remove the registration device with the transferbow from the mouth of the patient.


6.2.5 Reading off and documenting the angle values

Now, read the CE angle off the CE scale (9) and note the angle value on the corresponding data sheet or the order slip, respectively.

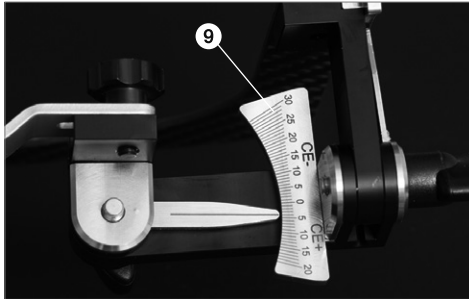
Now, also read the BP angle off the BP scale (6) and note the angle value on the corresponding data sheet or the order slip, respectively.

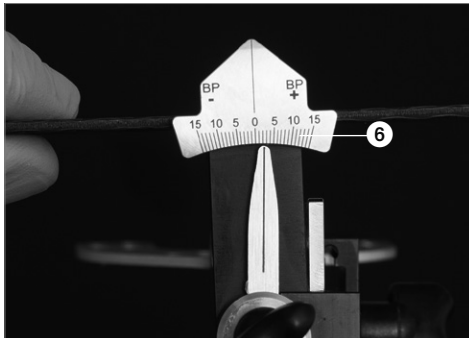6.2.6 Loosening the bite fork fixation screw

Loosen the bite fork fixation screw (10) and pull the bite fork (2) or centric tray off the bite fork adapter (13).

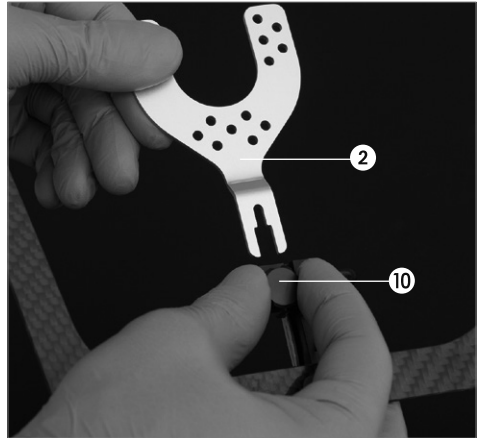6.2.7 Cleaning the components

Measurement is completed and the components can be cleaned and wipe disinfected. Please also observe Chapter 7.

7 Maintenance, Cleaning, Diagnosis

This chapter describes the user maintenance and cleaning procedures for UTS CAD. Only those tasks are listed that may be performed by dental professionals. All other tasks must be performed by qualified service personnel at a certified Ivoclar Vivadent Service Centre.

7.1 Monitoring and maintenance

The time for these maintenance procedures depends on the frequency of use and the working habits of the user. For that reason, the recommended times are only approximations.

7.2 Cleaning

Item	Frequency	Cleaning material
Bite fork that was in the mouth of the patient and came into contact with saliva or blood	After every use	Sterilize in the autoclave for 5.5 min at 134 °C using the pre-vacuum method. Observe the valid standards.

7.3 Maintenance

Avoid any contact with strong acids and solvents (e.g. MMA) to prevent damage of the surfaces.

Wipe disinfect after every use.

- Use only disinfectants suitable for stainless steel and light metals.
- Avoid any contact with strong acids, brines or solvents.
- Clean in an ultrasonic bath, with steam jets or in autoclaves.

8 What If...

This chapter will help you to recognize malfunctions and take appropriate measures or, if possible, to perform some minor repairs.

8.1 Technical malfunctions

Description	Instructions for Users	Measures
The surface of the device is scratched	Was the surface cleaned or manipulated with a sharp object?	Do not scratch the surface with a sharp object
The surface of the device is damaged or partially dissolved	Was the surface cleaned with acid or a solvent?	Do not use any acids or solvents for cleaning
The surface of the device is partially dissolved	Was the surface cleaned with an unsuitable disinfectant?	Do not use any unsuitable disinfectants
The joints for the angle setting can no longer be soundly locked	Have the components been correctly mounted after the device was disassembled?	Correctly assemble device according to the Operating Instructions or contact your Service Centre

8.2 Repair

Repairs may only be carried out by a certified Ivoclar Vivadent Service Centre. If repairs during the warranty period are not carried out by a certified Service Centre, the warranty will expire immediately. Please also refer to the corresponding warranty regulations.

9 Product Specifications

9.1 Delivery form

- 1 UTS CAD
- 1 Bite fork
- 1 Operating Instructions

Delivery forms may vary from country to country.

9.2 Technical data

Weight (approximately):

- Bite fork: 11 g
- UTS CAD basic bow: 20 g
- UTS CAD complete: 98 g

Dimensions:

- Basic bow side length: 200 mm
- Basic bow width, front: 155 mm, rear: 225 mm

Material:

- Metal components: Anodized aluminium, stainless steel
- Lateral arms: Carbon fibre compound, NBR (nitrile rubber)

Adjustment options:

- CE value range: +20° to -30°
- BP angle range: +/- 15°

10 Miscellaneous

Tips for transferring the angle values to the CAD software

Please observe the Instructions for Use for digital CAD dentures or the corresponding textbooks.

The apparatus has been developed solely for use in dentistry. Start-up and operation should be carried out strictly according to the Operating Instructions.

Liability cannot be accepted for damages resulting from misuse or failure to observe the Instructions. The user is solely responsible for testing the apparatus for its suitability for any purpose not explicitly stated in the Instructions. Descriptions and data constitute no warranty of attributes and are not binding.

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© Ivoclar Vivadent AG, Schaan/Liechtenstein
Ivoclar Vivadent Limited
Ivoclar Vivadent, Inc.