IvoBase[®] Hybrid

IvoBase[®] High Impact

Instructions for Use

Page 7

Gebrauchsinformation

Seite 11

Mode d'emploi

Page 15

Istruzioni d'uso

Pagina 19

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Página 23

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Página 27

Bruksanvisning

Sidan 31

Bruksanvisning

Side 35

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Sivu 39

Bruksanvisning

Side 43

Productinformatie

Pagina 47

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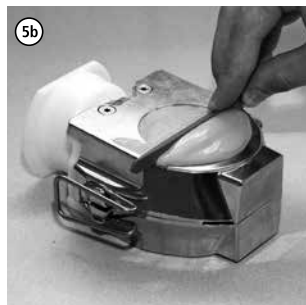
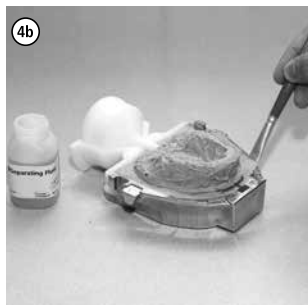
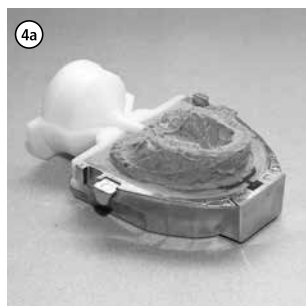
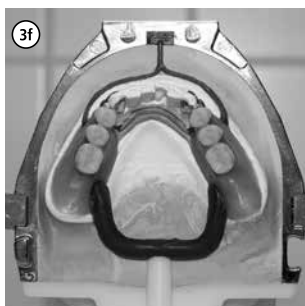
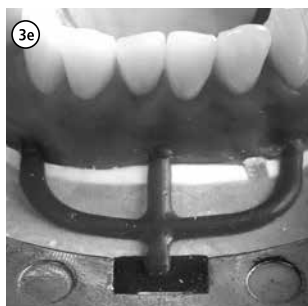
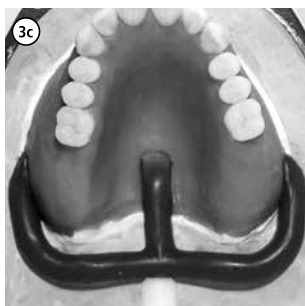
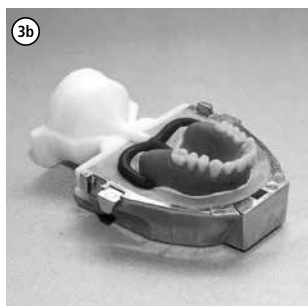
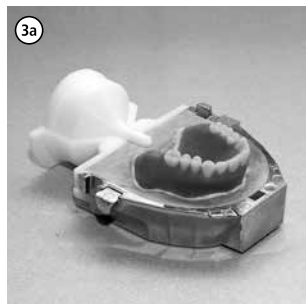
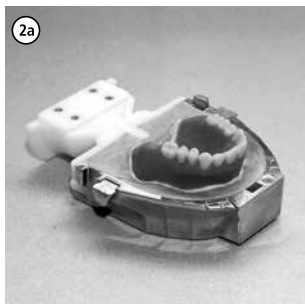
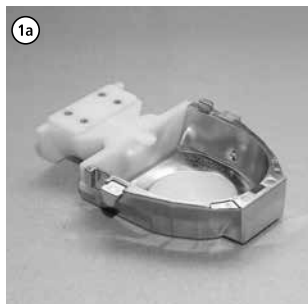
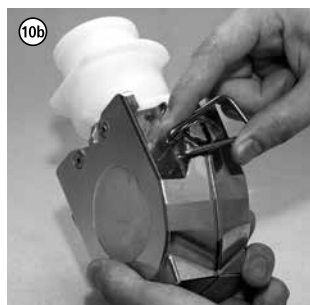
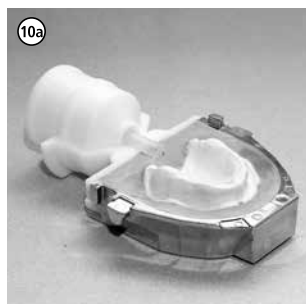
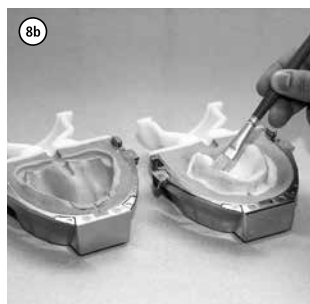
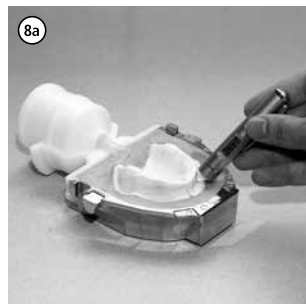
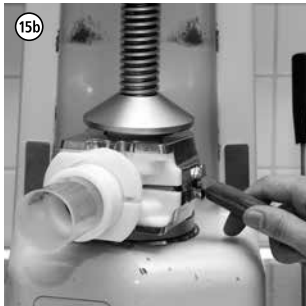
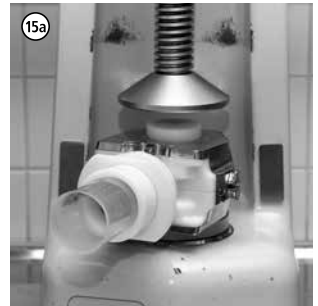
Instructions for Use

Please read these Instructions for Use carefully and familiarize yourself with the use of the material. Please also refer to the Operating Instructions of the Ivobase Injector.

Product description

The Ivobase system consists of a PMMA-based denture base material and the Ivobase Injector specifically developed for the material. The system enables fully automated and shrinkage-compensated polymerization of dentures.

The Ivobase material combines the advantages of heat- and auto-curing denture base materials. Even though it is subject to the requirements placed on auto-curing polymers, it also meets the stricter normative requirements for heat-curing polymers. Ivobase may thus be considered a hybrid material. In addition to the conventional hybrid material, Ivobase is also available as an impact-resistant material version: Ivobase High Impact.

The following shades are available:

- Ivobase Hybrid:
Pink, Pink-V, Preference, Clear, Pink-V Implant, Preference Implant
- Ivobase High Impact:
Pink, Pink-V, Preference, Pink-V Implant, Preference Implant, 34-V

Composition

- **Ivobase Hybrid**
Powder: Polymethyl methacrylate, citrate softener, initiator, pigments
Liquid: Methyl methacrylate, dimethacrylate, catalyst
- **Ivobase High Impact:**
Powder: High-impact-modified PMMA copolymer, PMMA copolymer, polymethyl methacrylate, citrate softener, initiator, pigments
Liquid: Methyl methacrylate, dimethacrylate, catalyst

Indication

- Complete denture prosthetics
- Partial denture prosthetics
- Combined denture prosthetics
- Hybrid denture prosthetics
- Implant prosthetics
- Relines
- Repairs
- Mouthguards

Contraindication

If the patient is known to be allergic to any ingredients. Intraoral use of unpolymerized material.

Side effects

In individual cases local allergic reactions to methacrylate materials have been reported.

Safety notes

- Monomer contains methyl methacrylate (MMA)
- MMA is highly flammable and irritating, flash point: +10 °C/50 °F.
- MMA is irritating to eyes, respiratory organs and skin.
- Skin contact may cause sensitization.
- Avoid skin contact with monomer and unpolymerized material. Commercial medical gloves do not provide protection against the sensitizing effect of methacrylates.
- Do not inhale vapours.
- Keep away from sources of ignition – do not smoke.
- Prevent monomer from reaching the sewage system.
- Take measures against electrostatic charge.

Notes on processing

- Use only Class III dental stone for model fabrication and investment.
- At warmer ambient temperatures, it is recommended to store the Ivobase wax components in the refrigerator at a temperature between 2 °C and 8 °C (36 °F–46 °F) for 5 minutes. This ensures that the wax components can be easily removed from the carrier foil.
- Thoroughly soak models in water for 5–10 minutes before investment.
- Thoroughly isolate inner flask surfaces before investment.
- Clean stone surfaces and isolate them with a thin layer of Ivoclar Vivadent Separating Fluid.
- Closely observe polymerization and cooling times.
- Contact of solvents or monomer with polymerized denture base material may lead to white discolouration.
- When using cervical isolations, Ivoclar Vivadent recommends A-silicones with a hardness of >65 Shore A (e.g. Flexistone Plus/Virtual® Heavy Body Fast). Condensation cross-linking silicones (C-silicones) may have a damaging effect on the denture base material.
- The products may only be used and processed by trained personnel.
- Post-injection of a second capsule of Ivobase or SR Ivocap® material into an already filled flask is not permissible.

- The remaining material of an already injected IvoBase capsule cannot be used again.
- Dried stone must be soaked in water prior to injection.
- A contaminated flask may lead to smoke development when heated.
- Store in a humid environment from the completion of the dentures to their incorporation.

Processing

Preparing for injection

1. Preparing the flask

- ①a Isolate the inner surfaces of the clean flask halves with a thin layer of petroleum jelly.

Place the flask lid, access former half and the filter wax component (B1) in one of the flask halves from the inside.

2. Investing the model

- ②a Invest the well-soaked and isolated model in one of the flask halves using conventional dental stone (Class III). Isolate the side and bottom surfaces of the model exclusively with Ivoclar Vivadent Separating Fluid (do not use petroleum jelly or soap). Skim off the excess stone in such a way that it is flush with the model and the edge of the flask housing.

Important: The stone edge must be flush with the access former half.

3. Positioning the wax components

The wax of the wax components (B) can be secured on the contouring and the model by means of wax. Make sure that no stone can flow below the wax components.

- ③a Allow the stone to set thoroughly. Remove the access former half and place the access former full.

– Place the injection wax component (B2)

- ③b Press the injection wax component (B2) onto the incline of the access former full.

Maxillary complete dentures

- ③c Position a sprue at the centre of the palate roof and on the vestibular side of each upper cusp and press the wax onto the model.

Mandibular complete dentures

- ③d Place one sprue each to the lingual area of the retromolar triangle. Cut off the centre sprue of the injection wax component (B2).

Partial dentures

Place a sprue at the dorsal end of each denture saddle.

– Place the aeration wax component (B3)

- ③e Press the aeration wax component with the aeration channel to the centre of the filter wax component.

Important: The aeration channels must not come into contact with the flask housing when being pressed into place in order to ensure the tightness of the flask.

Mandibular and maxillary complete dentures:

- ③f Press the three channels of the aeration wax component (B3) to the anterior area of the contouring.

Partial dentures:

Press the three channels of the aeration wax component (B3) to the anterior stops of the denture saddles. If the aeration wax component (B3) is not long enough, the connection between the denture saddle and the filter wax component has to be completed with a wax sprue observing the shortest possible distance.

4. Investment with silicone

- ④a Cover the vestibular and palatal and/or lingual surfaces of the contoured dentures with silicone. Do not cover the occlusal surfaces and incisal edges of the teeth. Apply a retentive pattern in the surface of the silicone before it sets or use retention sand to secure the silicone in the stone.

Note: For optimum results, the use of a flowable, addition cross-linking silicone (A-silicone with a Shore-A hardness of >65) is recommended (e.g. Flexistone Plus or Virtual Heavy Body Fast).

- ④b Isolate stone-to-stone contact surfaces with a layer of Separating Fluid (F) (do not use petroleum jelly or soap).

5. Casting the counter model

- ⑤a Position the upper flask half and lock using the locking clasps.

- ⑤b Mix conventional dental stone (Class III) according to the instructions of the manufacturer and have it flow in flowable consistency into the opening of the flask without any air bubbles with the help of a shaker. Skim off excess stone with the IvoBase spatula so that no stone protrudes from the flask opening.

6. Boiling out the model

- 6a) Once the stone has entirely set, heat the flask in a water bath at approximately 90 °C (194 °F) for approximately 5–8 minutes. After that, open the flask, remove the access former full and roughly remove the wax using a wax knife. Now thoroughly boil out the inner surfaces with clean, boiling water.

Note: Use only pure water without additives (e.g. wax solvent, cleaning agent). Additives may result in a discolouration of the dentures.

7. Teeth

- 7a) For an optimum bond, tooth necks and the basal tooth surfaces must be absolutely free of wax. Slightly blast the basal and cervical tooth surfaces (100 µm Al₂O₃, 2 bar / 29 psi) or roughen them with a dental tungsten carbide bur.

8. Isolating the stone surfaces

- 8a) **Important:** Before isolation, cool the flask halves to room temperature (<30 °C / 86 °F) using air or clean, cold water. The exact flask temperature can be checked with the help of an IR thermometer.
- 8b) Isolate the clean surfaces of the plaster model with Separating Fluid (F). Apply a second coat after 5 minutes. Pooling (e.g. in the gingivo-buccal fold) must be prevented.

9. Placing the aeration filter

- 9a) Insert the aeration filter (E) flush into the corresponding recess in the flask.

10. Placing funnel and centring insert

- 10a) Push the IvoBase funnel (D) into the centring insert until it snaps into place so that it is locked. Place the centring insert with the inserted funnel (D) in the lower flask half. Make sure that the isolation has dried.
- 10b) Assemble the flask halves, check for contamination and lock the halves using the locking clasps. Narrow gaps are closed by the pressure of the locking clasps.

Injection and polymerization

The operation of the IvoBase Injector is described in detail in the corresponding Operating Instructions. The IvoBase Injector must be brought to its stand-by state before the material is mixed.

11. Mixing the material

One capsule is sufficient for the injection of a maxillary or mandibular complete denture.

Contents of the capsule:	Monomer	Polymer
IvoBase Hybrid	20 ml	34 g
IvoBase High Impact	20 ml	30 g

With very large denture bodies, for which more resin is required, mix the material of two capsules in a separate mixing cup and fill the required quantity back into one capsule. Subsequently, injection and polymerization are conducted as described below.

- 11a) Remove the monomer container (A4) from the capsule. Remove the capsule seal (A1).
- 11b) Twist off the sealing cap of the monomer container (A4) and pour the liquid into the open capsule (A2).
- 11c) Use the spatula (C) to mix the monomer and polymer in the capsule to a homogeneous mixture (for approximately 20–30 seconds).

12. Loading the Injector

Proceed as follows within 5 minutes after mixing:

- 12a) Place the open capsule (A) on an even surface and press the flask onto the capsule (A).

Note: Do not tilt the flask with the capsule prior to polymerization to prevent the material from flowing out of the capsule (A). When the flask is lifted, the capsule (A) must be located securely in the centering insert.

- 12b) Open the door of the Injector. Slide the flask with the capsule (A) into the polymerization chamber until it stops via the flask holder. Close the door of the Injector.

13. Starting the program

- 13a) Select the corresponding program (see Operating Instructions of the Injector for a program overview).
- 13b) If desired, activate the RMR key in order to reduce the residual monomer content to <1%.
- 13c) Press "Start". The fully automated injection process begins. The remaining process time is indicated in the display.
- 13d) At the end of the program, an acoustic signal sounds. Use the thermal glove to remove the hot flask from the Injector by holding it by the isolation shoulder. The Injector is ready for the next polymerization immediately after the the flask is removed.
- Note:** The flask may still be very hot when removed from the Injector. Burn hazard!

14. Cooling the flask

Cool the flask under cold running water (water temperature <15 °C / 59 °F) for at least 15–25 minutes. As an alternative, cooling may also take place in the Injector (e.g. over night).

Important: If the flask is not cooled down and opened at too high a temperature, inaccuracies of fit may result.

15. Divesting the dentures

Remove the locking clasps.

- 15a) Place the flask under a dental press. Secure the divesting aid between stone and pressing plate.
- 15b) Slide a broad screw driver or plaster knife into the joint between the two flask halves and lever the upper flask half upwards. Remove the flask from the press, turn it around and repeat the procedure.
- 15c) Remove the stone core with the capsule and separate the capsule at the sprue using a separating disk bur or saw. Divest the denture from the stone core as usual with stone tongs.

16. Finishing the dentures

The polymerized IvoBase material is finished with cross-cut tungsten carbide burs, sandpaper and polishing media as it is usual for denture base materials.

Repair / relining

IvoBase can be repaired with itself as well as with ProBase® Cold auto-curing polymer.

Repair / extension with IvoBase Hybrid / High Impact

Basically, an IvoBase denture can be repaired without investment in the IvoBase flask. The basic procedure corresponds to that used for repairing conventional auto-curing denture base materials.

Observe the mixing ratio as follows:

	Polymer	Monomer
IvoBase Hybrid	9 g	5 g
IvoBase High Impact	8 g	5 g

The polymer and monomer must be homogeneously mixed (for approximately 20 seconds).

Please observe the following processing times (at 23 °C / 73 °F) room temperature):

- Duration of the dough phase: 30 seconds
- Duration of the pouring phase: 1–2 minutes
- Duration of the modelling phase: 3–5 minutes

The maximum total processing time after mixing is 10 minutes.

Polymerize the material in a pressure pot for 20 minutes at 55 °C / 131 °F at 2.5 bar / 36 psi pressure.

Relining with IvoBase Hybrid / High Impact

The relining of an IvoBase denture may be carried out with or without investment in the IvoBase flask.

- a) Relining without investment in the IvoBase flask
The procedure without investment in the IvoBase flask is the same as that for relining with conventional auto-curing denture base materials.
- b) Relining with investment in the IvoBase flask
Investment and completion for the relining of a denture is basically the same as for the completion of a waxed-up denture. Make sure that
 - the sprues and aeration channels are attached to the impression of the relining.
 - the area to be relined demonstrates a thickness of at least 2 mm.

Repair / extension relining with ProBase Cold

The procedure corresponds to that used for repairing conventional auto-curing denture base materials. Please refer to the ProBase Cold Instructions for Use for details.

Storage instructions

- Store material in a cool, dark and well ventilated place.
- Storage temperature: 2–28 °C / 36–82 °F.
- Do not use the products after the indicated expiration date.
- Keep material out of children's reach.

Classification according to EN ISO 20795-1, Type 2, Class 1

The product meets the requirements of EN ISO 20795-1:2013.

This material has been developed solely for use in dentistry. Processing should be carried out strictly according to the Instructions for Use. Liability cannot be accepted for damages resulting from failure to observe the Instructions or the stipulated area of application. The user is responsible for testing the material for its suitability and use for any purpose not explicitly stated in the Instructions. Descriptions and data constitute no warranty of attributes.

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Ivoclar Vivadent Limited
Ivoclar Vivadent, Inc.