



Tetric® PowerFill | Tetric® PowerFlow

1

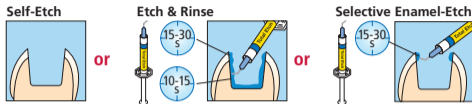


2

Adhese® Universal
VivaPen®

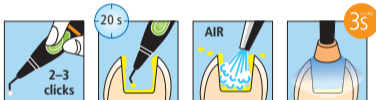

3

Tetric®
PowerFlow

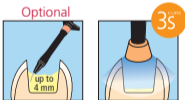

4

Tetric®
PowerFill



Bluephase®
PowerCure



	Tetric PowerFill / Tetric PowerFlow	Adhese Universal
500–900 mW/cm ²	20 s	10 s
900–1 400 mW/cm ²	10 s	10 s
1 800–2 200 mW/cm ²	5 s	–
2 700–3 300 mW/cm ²	3 s	3 s

Direct posterior restorations, light-cured from occlusal aspect
Direkte Seitenzahnrestorationen, okklusale Belichtung



Date information prepared:
2018-07-10
693240/Rev. 0

CE 0123

Ivoclar Vivadent AG
FL-9494 Schaan/Liechtenstein

ivoclar
vivadent
clinical